



TRAVELER 4you Ergo

PERFECT FIT
MEETS THE SMALLEST PACK SIZE.

The TRAVELER 4you Ergo is your individual active folding wheelchair with a fine frame geometry and ergonomic seat fold. You can choose 12 frame parameters, such as knee angle, frame retraction, wheelbase, or the position of the caster wheel bearing block. The frame is made according to your wishes as a one-off. With a compact pack size of 17 x 35 x 58 cm, you can even take your TRAVELER 4you Ergo with you in a sports bag. Many settings, such as tipping point, seat height, seat angle or back angle, can also be adjusted to your needs afterwards.




#Reliability

Thanks to its excellent stability, you can combine the TRAVELER 4you Ergo with adaptive bikes, traction devices and steering attachments as well as other adaptations.



#Individuality

With 12 selectable frame parameters, each frame is unique. The frame stem can include a frame retraction.

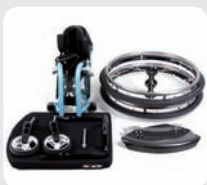

#Independence

With its small pack size, the folding backrest and tool-free detachable assemblies, the TRAVELER 4you Ergo is easy to transport.



#MoreMobility

By simply folding it, you can reduce the seat width (sitting in a wheelchair) by 15% or 5-10 cm to be able to pass through narrow passages effortlessly.



INDIVIDUALLY ADAPTED, ERGONOMICALLY DESIGNED.

Your **TRAVELER 4you Ergo** is manufactured individually according to your needs, from the frame dimensions to the equipment. You can choose from options such as footrests, brakes, push handles, seat and back systems, clothing guards, armrests, and drive and caster wheels to tailor your TRAVELER perfectly to your needs.

UNIQUE, SELECTABLE OPTIONS FROM PRO ACTIV:



The 7-way angle-adjustable and foldable backrest offers maximum flexibility for the ideal sitting position.



The footrest, which can be folded towards the back with a spring locking mechanism, is easy to handle and flexible in everyday life.



EQUIPMENT PACKAGE "PURE":

With this package you can achieve a particularly light total weight from 10 kg, which sets new standards in terms of lightness. It offers innovative, weight-optimised features such as carbon side panels with integrated clothing protection, a "light" adjustable back with back padding and a lightweight footrest with rounded footplate support without protrusion.




TECHNICAL SPECIFICATIONS

Seat width:	33 - 52 cm (in 1 cm steps) 33 - 42 cm with "pure" equipment
Seat depth:	35 - 48 cm (in 1 cm steps)
Back height:	22 - 48 cm (in 1 cm steps)
Lower leg length:	37 - 51 cm (in 1 cm steps)
Rear seat height:	39 - 47 cm (in 1 cm steps)
Front seat height:	42 - 52 cm (in 1 cm steps)
Tipping point:	3 - 11 cm (in 1 cm steps)
Total product weight:	Starting at 11.2 kg Starting at 10.0 kg with "pure" equipment
Max. standard load weight:	Up to 120 kg

Data subject to design changes.

FOLLOW US ON SOCIAL MEDIA

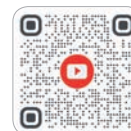

PROACTIV[®]

PRO ACTIV Reha-Technik GmbH

Im Hofstätt 11
72359 Dotternhausen
info@proactiv-gmbh.de

www.proactiv-gmbh.com

Rev. 09/2024

Subject to technical modifications, printing errors & colour deviations.