



SPEEDY F1

PRECISION AND SPEED FOR EVERY GAME.

With 15 selectable frame parameters such as knee angle, frame retraction, wheelbase and lower leg width at the top and bottom, the **SPEEDY F1** is the individualist in its class. Thanks to the robust frame geometry and the high-quality 40 mm oversized tubes, you get excellent riding characteristics and manoeuvrability. The perfect fit and weight distribution ensure that your muscle power is optimally converted into movement - for precise and demanding manoeuvres.



SPEEDY F1 4basket

Experience easy turning, fast responsiveness and acceleration in a modern design. Various front bars for offence and defence and rear ram bars and anti-tipping supports adapt to your playing style and playing position.



SPEEDY F1 4tennis

The tennis frame offers you optimal playing clearance with torsional rigidity and direct response. The height-adjustable anti-tipping support can be quickly adapted to the ground conditions and your playing habits.



SPEEDY F1 4badminton


The badminton frame ensures aggressive positioning and stability with widely protruding tilt supports and a wide wheelbase. In combination with a large camber, these features results in extreme manoeuvrability and optimal acceleration behaviour.



YOUR SPEEDY F1 FOR TOP ATHLETIC PERFORMANCE.

Your **SPEEDY F1** is individually manufactured according to your needs in terms of frame dimensions and equipment. You can choose from a wide range of options, including specific details such as a high handrim attachment for optimal power transmission and high acceleration.

UNIQUE, SELECTABLE OPTIONS FROM PRO ACTIV:

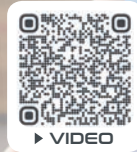

For users with hemiplegia or limited mobility on one side of the body, there is the manoeuvrable one-hand drive with double handrim up to 12° camber.



Choose between anti-tipping supports on both sides or a centre anti-tipping support for optimal stability.



As an option, different footplate supports such as “front edge” or “box shape” can be selected, which provide secure positioning for the feet.

TECHNICAL SPECIFICATIONS

Seat width:	Custom
Seat depth:	Custom
Back height:	Custom
Lower leg length:	Custom
Rear seat height:	Custom
Front seat height:	Custom
Tipping point:	Custom
Total product weight:	Starting at 7.8 kg
Max. standard load weight:	Up to 120 kg

FOLLOW US ON SOCIAL MEDIA



PROACTIV[®]

PRO ACTIV Reha-Technik GmbH

Im Hofstätt 11
72359 Dotternhausen
info@proactiv-gmbh.de

www.proactiv-gmbh.com